


HORIZON™ 627 LUMBAR, 631 LSO LoPro, 637 LSO

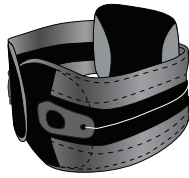Horizon™ 627 Lumbar

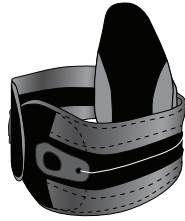Horizon™ 631 LSO Lo Pro

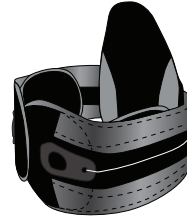
Horizon™ 637 LSO

Distributed By

Zynex Customer Service

10000 Park Meadows Drive
Lone Tree, CO 8024
P: 1-(800) 495-6670
E: customersupport@zynex.com

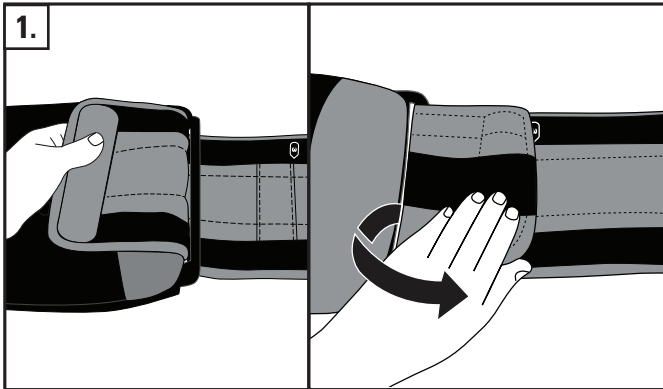1.

**ENGLISH
SIZE ADJUSTMENT**

To change the size of the Horizon, peel up the inside end of the side panel on each side. Adjust to proper fit and reattach.

The Horizon™ Brace sizing follows closely to standard pants sizing. Size ranges are only guidelines. Your experience and personal assessment could make the difference in choosing the correct size.

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
	1	2	3	4
HORIZON 627 LUMBAR	24 to 34 in. (66 to 86 cm)	32 to 40 in. (81 to 102 cm)	38 to 46 in. (97 to 117 cm)	44 to 50 in. (112 to 127 cm)
Men's Pant Size	26 - 34	30 - 38	34 - 42	38 - 46
Women's Pant Size	0 - 4	6 - 10	12 - 16	18 - 24

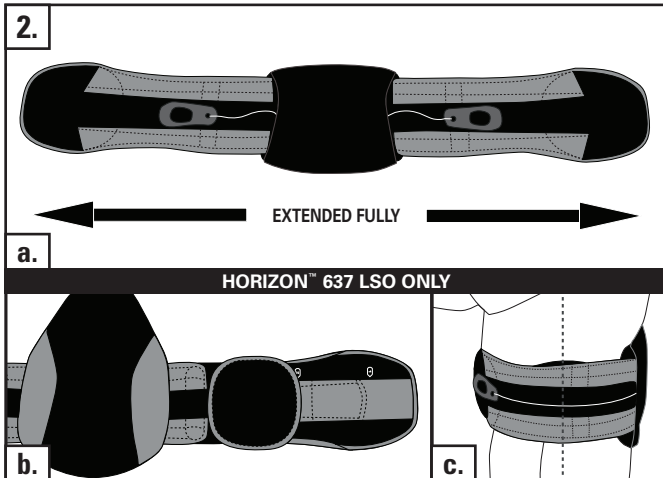2.

a. EXTENDED FULLY

HORIZON™ 637 LSO ONLY

b.

c.

**ENGLISH
PREPARATION**

a. Place the pull-tabs 2 inches (5 cm) from the plastic in the back of the Horizon™ brace. Pull the ends of the Horizon™ brace away from each other ensuring that it is fully extended.

HORIZON™ 637 LSO ONLY
b. With brace extended fully, place one lateral panel on midline on each side panel of the brace. Press down to secure to brace.

c. When Horizon™ brace is properly fit to the patient, the lateral panels should fit at patient's midline.

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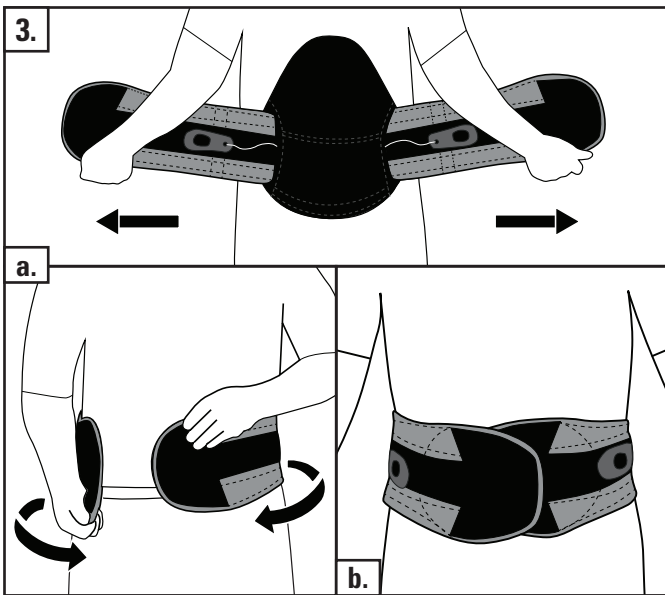
ENGLISH
APPLICATION

a. Apply the fully extended brace making sure the back panel is centered and low on the back. Apply the brace tightly around the abdomen, the right side overlapping the left.

b. When applied properly, both the upper and lower edge of the brace should be snug. If there is any gap, reapply the brace and angle the ends to insure a snug fit at both the upper and lower edge.

Note: For certain patients, the ends must be steeply angled.

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
LIMITED WARRANTY
Aspen Medical Products, Irvine, CA 92618, warrants to the user who originally purchases this product that it is free from defects in material and workmanship. The sole obligation of Aspen Medical Products in the event of breach of warranty shall be to repair or replace the defective product or part(s).
Aspen Medical Products shall have no obligation under this limited warranty in the event:
(a) The product was not purchased from Aspen Medical Products or through its authorized channels of distribution;
(b) The product is altered;
(c) Any parts not supplied by Aspen Medical Products are inserted into the product; or
(d) The product is not used in accordance with the Aspen Medical Products, Instructions for Use.

THE FOREGOING IS THE SOLE AND EXCLUSIVE REMEDY FOR FAILURE IN SERVICE OF, OR DEFECT IN, THE PRODUCT. ASPEN MEDICAL PRODUCTS SHALL NOT BE LIABLE UNDER THIS OR ANY IMPLIED WARRANTY FOR ANY DIRECT SPECIAL, INCIDENTAL, OR CONSEQUENTIAL DAMAGES. THIS WARRANTY IS IN LIEU OF ALL OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING THE WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE OR USE, AND ALL OBLIGATIONS OR LIABILITIES ON THE PART OF ASPEN MEDICAL PRODUCTS FOR DAMAGES ARISING OUT OF OR IN CONNECTION WITH THE USE OF THE PRODUCTS, WHICH ARE HEREBY DISCLAIMED AND EXCLUDED BY ASPEN MEDICAL PRODUCTS.
This warranty gives you specific legal rights and you may have other rights which vary from state to state.
Warning: To ensure proper use, please review all material.
A listing of authorized distributors in your area is available upon request.

READ INSTRUCTIONS BEFORE USE. Proper training in the use of this device should take place before it is applied. These directions are guidelines only and are not offered as medical recommendations.


THIS PRODUCT IS INTENDED FOR APPLICATION BY LICENSED HEALTH CARE PRACTITIONERS AS DIRECTED BY A PHYSICIAN OR OTHER QUALIFIED MEDICAL AUTHORITY. PERMANENT INJURY MAY RESULT FROM USE OTHER THAN AS DIRECTED BY A PHYSICIAN OR OTHER MEDICAL AUTHORITY. ASPEN MEDICAL PRODUCTS SPECIFICALLY DISCLAIMS ANY OBLIGATION OR LIABILITY FOR INJURY OR DAMAGE BY REASON OF ANY APPLICATION OF THIS PRODUCT. PLEASE SEE THE ACCOMPANYING LIMITED WARRANTY FOR FURTHER INFORMATION IN THIS REGARD.

HORIZON™ LATERAL EXTENSION PANELS

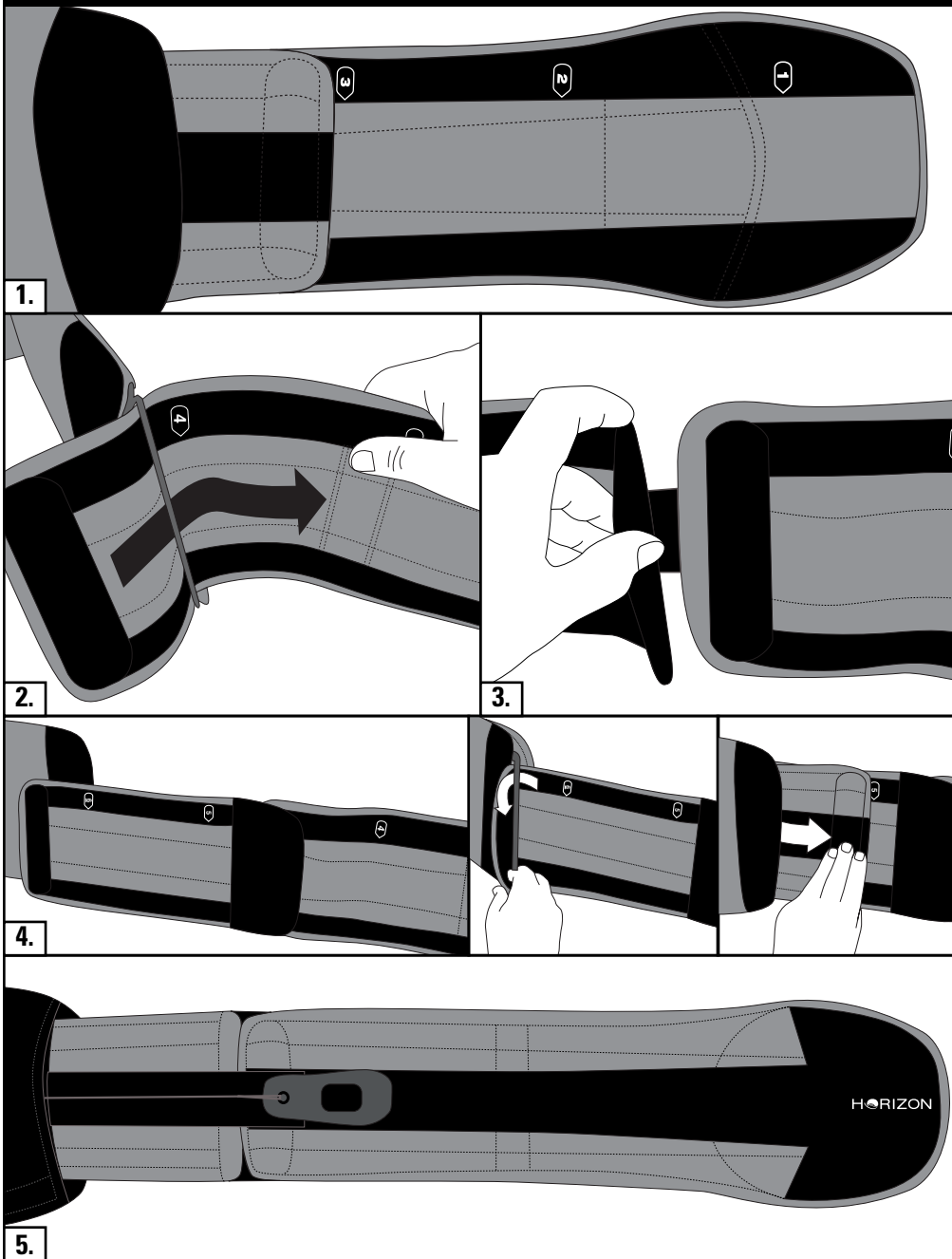HORIZON™

LATERAL EXTENSION PANELS

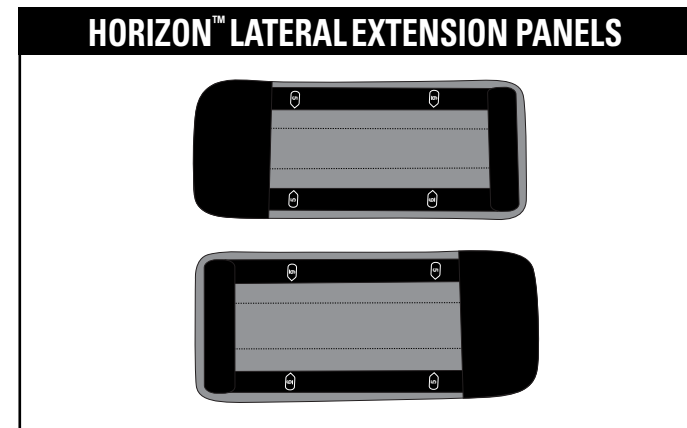	3	4	5	6
EXTENSION ACCESSORY	50 to 58 in. (127 to 145 cm)	56 to 62 in. (142 to 160 cm)	60 to 68 in. (157 to 175 cm)	66 to 70 in. (173 to 193 cm)
Men's Pant Size	42 - 50	46 - 54	50 - 58	54 - 62
Women's Pant Size	26 - 32	34 - 40	42 - 48	50 - 56

Aspen Medical Products

6481 Oak Canyon, Irvine, CA 92618-5202, USA

800-295-2776 • 949-681-0200

aspennmp.com