


E-Wave™

NeuroMuscular Electrical Stimulator

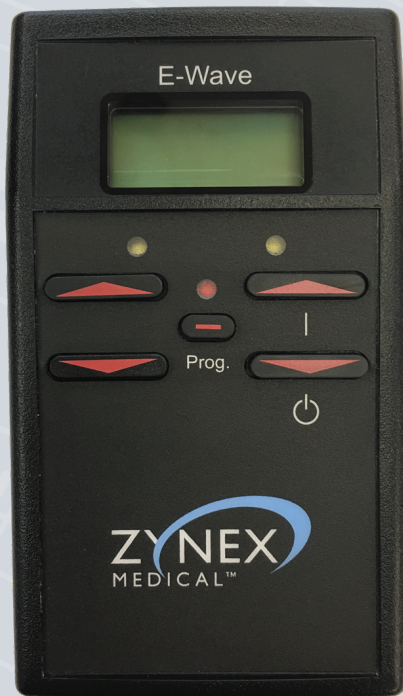

Rx Only - Made in USA

New Technology

The E-Wave™ is an easy-to-use NeuroMuscular Electrical Stimulator (NMES) with unique features designed for high performance and a wide variety of options.

The E-Wave™ is one of the safest stimulators in the world with safety features such as electrode alarm, level locking and a guaranteed constant current.

The E-Wave uses modern state-of-the-art technology to accurately deliver the desired intensity and frequency of electrical current. The E-Wave™ is easy to program with a digital on-screen display, with step-by-step options which guide the user through the set-up of muscle stimulation protocols.

- Fully programmable, small battery powered device
- Adjustable treatment timer and compliance meter
- Optional external trigger input to be used with a heel or hand switch

Adjustable Parameters

| | |
|--------------|--|
| Amplitude: | 0-100 mA |
| Frequency: | 4-100 Hz |
| Pulse-Width: | 50-300 μ sec |
| Ramp Up: | 0.1 to 6 Sec. |
| On-Time: | 0.5 to 30 Sec. |
| Ramp Down: | 0.1 to 6 Sec. |
| Off-Time: | 0.1 to 60 Sec. |
| TX Timer: | Continuous or 10-100 Minutes |
| Waveforms: | Symmetrical Biphasic (AC) or Monophasic (DC) |



External Hand Switch



External Heel Switch

Compatible Accessories

The remote switch allows either the patient or therapist the ability to trigger a muscle contraction manually.



Form 30087 Rev. 4